



Mirau Interferometer

(10X, 20X, 50X)

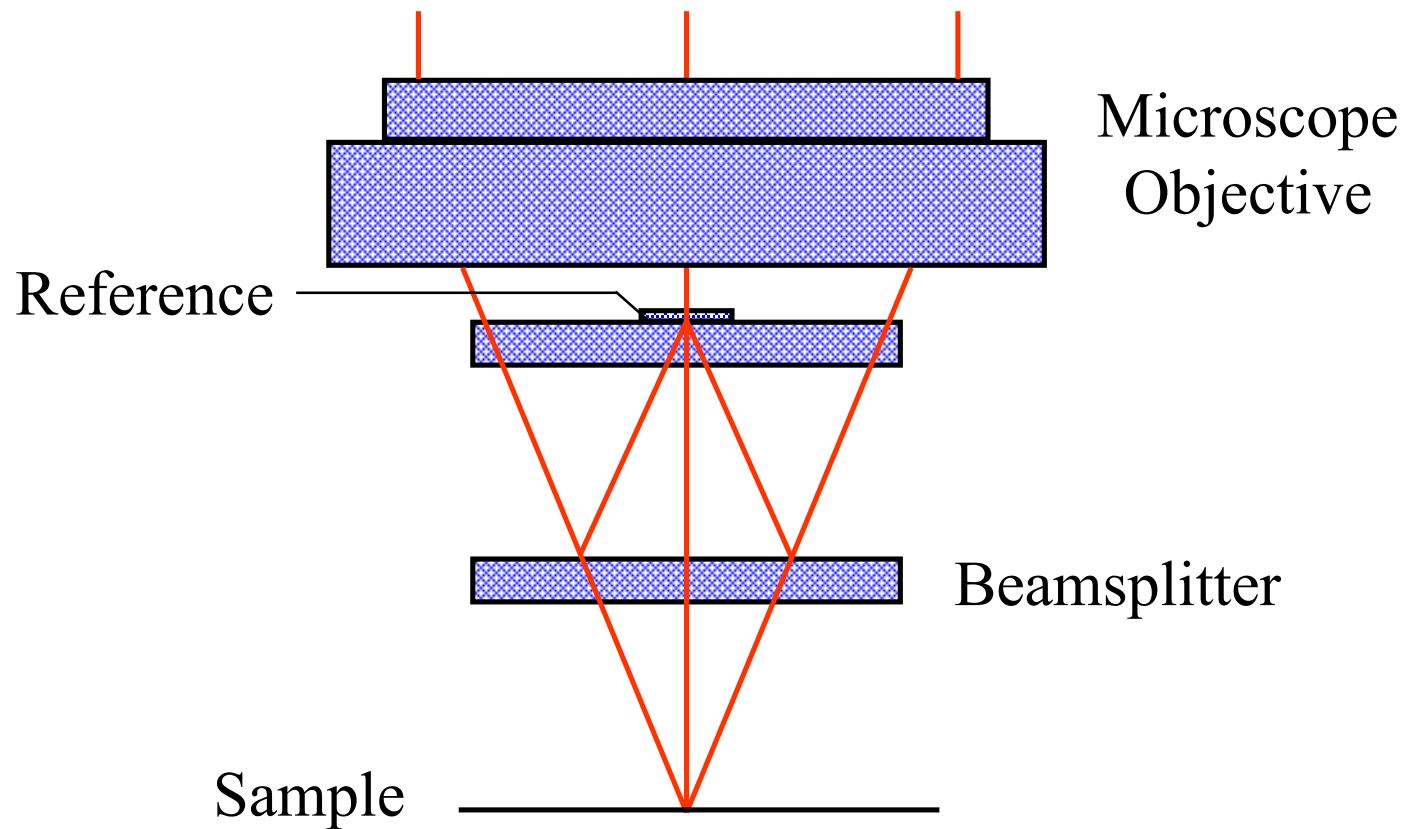


Michelson Interferometer

(1.5X, 2.5X, 5X)

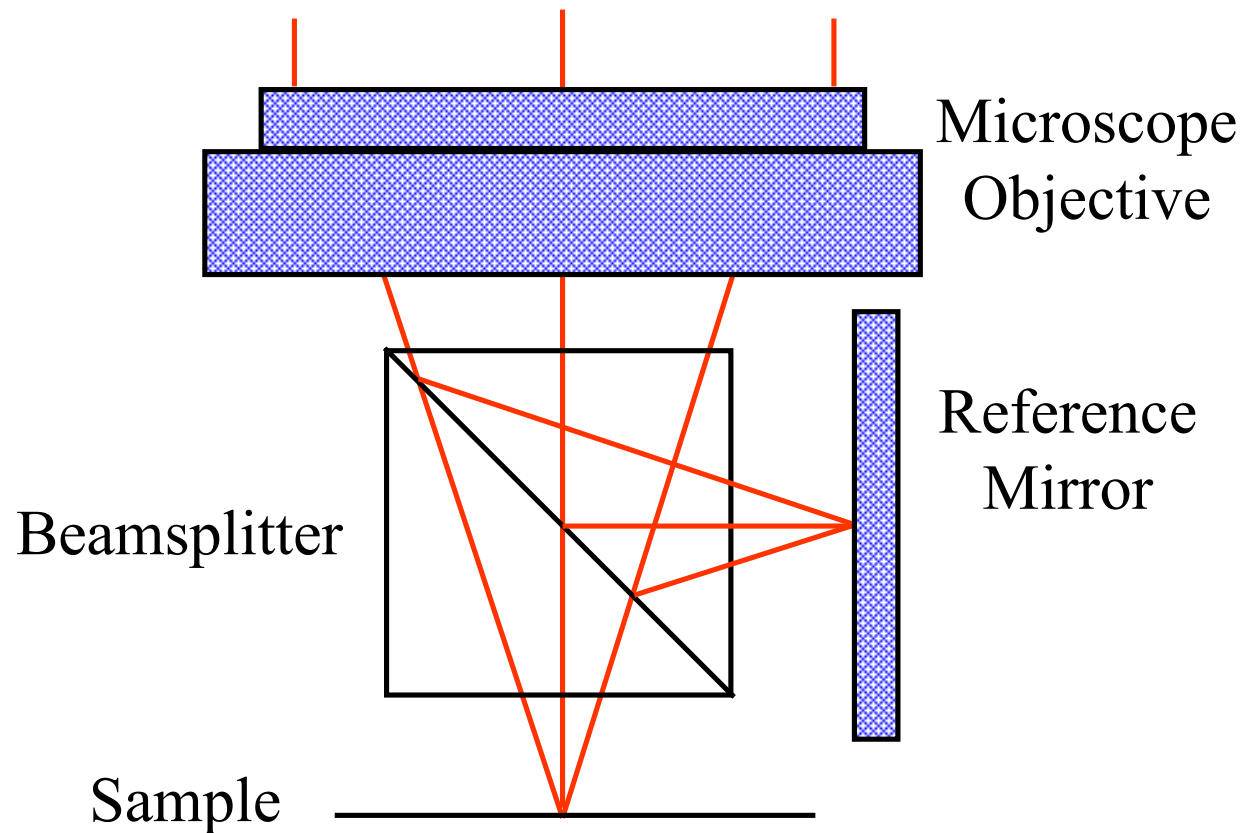

Linnik Interferometer

(100 X, NA 0.95)

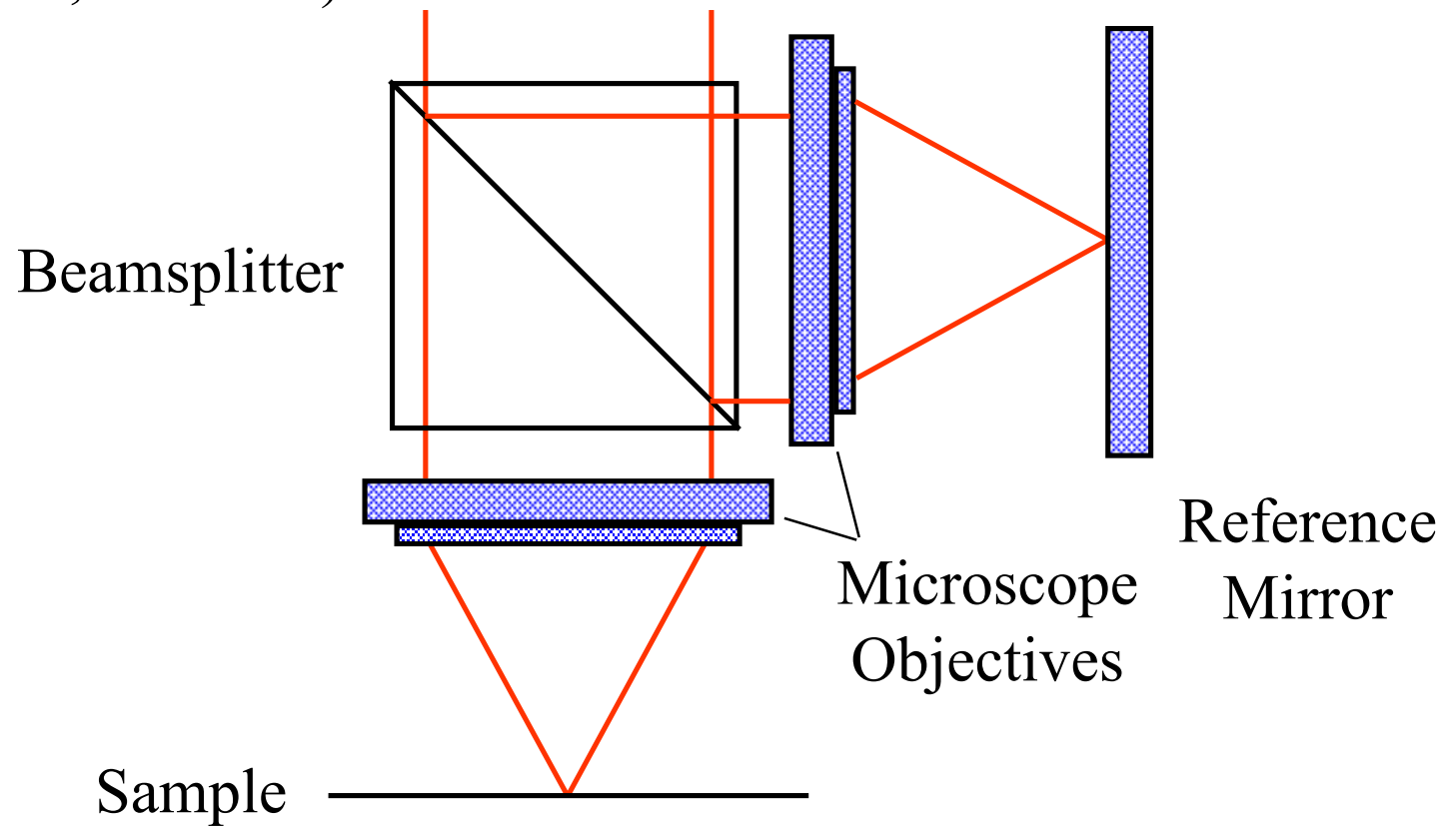

Interference Objectives

- **Mirau**
 - Medium magnification
 - Central obscuration
 - Limited numerical aperture
- **Michelson**
 - Low magnification, large field-of-view
 - Beamsplitter limits working distance
 - No central obscuration
- **Linnik**
 - Large numerical aperture, large magnification
 - Beamsplitter does not limit working distance
 - Expensive, matched objectives